

amplivox

Model 116

Operation Manual



ABOUT THIS MANUAL

READ THIS OPERATING MANUAL BEFORE ATTEMPTING TO USE THE INSTRUMENT.

This manual is valid for the Model 116 (applies from firmware version 4v47 onwards – please refer to section 1.3).

Amplivox Ltd.
10393 West 70th Street
Eden Prairie
MN 55344

This product is manufactured by:

Amplivox Ltd.
3800 Parkside, Solihull Parkway,
Birmingham Business Park,
Birmingham, West Midlands,
B37 7YG

www.amplivox.com

For all enquiries please contact us under:

Tel: 888 941 4208

Tel: +44 (0)1865 880846

Fax: 952 903 4100

Fax: +44 (0)1865 880426

info@amplivox.us

hello@amplivox.com

CE 0088

For supply in US only

Caution: Federal Law restricts this device to sale by or on the order of a licenced medical professional.

Contents

| | |
|---|-----------|
| 1 Introduction | 1 |
| 1.1 Intended applications | 1 |
| 1.2 Unpacking | 1 |
| 1.3 Firmware version | 1 |
| 1.4 Standard contents | 2 |
| 1.5 Optional accessories | 2 |
| 2 Important Safety Instructions | 2 |
| 2.1 Precautions | 2 |
| 2.2 Electromagnetic compatibility (EMC) considerations | 3 |
| 2.3 Power supply options | 3 |
| 2.4 Audiometer connections | 5 |
| 2.5 Data transfer to a printer | 6 |
| 2.6 Data transfer to a computer | 6 |
| 3 Using the Audiometer | 7 |
| 3.1 Switching the audiometer on and off | 7 |
| 3.2 Testing the patient response switch | 7 |
| 3.3 Audiometer display | 7 |
| 3.4 Audiometer controls | 7 |
| 3.5 Threshold Retention Function | 10 |
| 3.6 Printing audiograms | 11 |
| 3.7 Data transfer to Audibase or ampliSuite | 11 |
| 4 Suggested Sequence of Operation and Test Procedure | 12 |
| 5 Specification | 15 |
| 5.1 Output data | 15 |
| 5.2 Maximum hearing levels provided at each frequency | 15 |
| 5.3 Physical Data | 15 |
| 5.4 Equipment classification | 16 |
| 6 Symbols | 17 |
| 7 Technical Information | 18 |
| 8 Routine Maintenance | 19 |
| 8.1 Audiometer maintenance | 19 |
| 8.2 Transducer maintenance | 19 |
| 8.3 Mains adapter maintenance | 20 |
| 8.4 Batteries | 20 |
| 9 Instrument Storage and Transportation | 20 |
| 10 Calibration and Repair of the Instrument | 21 |
| 11 Guarantee | 21 |
| 12 Ordering Consumables and Accessories | 22 |
| 13 Disposal Information | 23 |
| Appendix 1 - EMC Guidance & Manufacturer's Declaration | 24 |
| Appendix 2 - Use with Non-medical Electrical Equipment | 31 |

1 Introduction

Thank you for purchasing an Amplivox audiometer. The Amplivox Model 116 is a manual screening audiometer that will give many years of reliable service if treated with care.

1.1 Intended applications

The Model 116 screening audiometer is designed for use by general practitioners, occupational health staff and child health professionals and is the ideal instrument for primary care groups, schools and industry. The audiometer is not intended for use by hearing specialists to determine the full extent and scope of a patient's hearing deficiency.

The instrument is completely portable, and if required may be specified to operate from integral internal batteries (see Section 2.3 below). Test results may be printed using the specified printer option or transferred to a PC running the Amplivox Audibase or ampliSuite applications.

1.2 Unpacking

Open the shipping carton and carefully remove all the equipment. Check against the delivery note that all the accessories ordered have been included with your audiometer. If anything is missing, please contact Amplivox Customer Support (+44 1865 880846; sales@amplivox.ltd.uk). If you have purchased from a distributor you should contact them directly.

Please retain the shipping carton and packing materials as the audiometer will need calibrating on an annual basis and should be returned to Amplivox in its original shipping carton.

1.3 Firmware version

This operating manual is for firmware versions 4v47 onwards. To check the version of firmware on your audiometer press and hold the MENU button followed by the TALKOVER button.



1.4 Standard contents

| | |
|-------------------------------|-------------------------|
| Model 116 Audiometer | Audiometric headset |
| Carrying case | Patient response switch |
| Mains adapter | Audiogram cards |
| Operating manual & ampliSuite | Calibration certificate |

1.5 Optional accessories

| | |
|--|----------------------------|
| Battery power function | Additional audiogram cards |
| Audibase software | USB Cable |
| Printer(s) | Printer cable(s) |
| Audiocups (noise reducing earphone enclosures) | |

2 Important Safety Instructions



The Model 116 instrument must be used only by practitioners qualified to perform audiometric tests. It is intended for use as a screening tool.

2.1 Precautions

READ THIS OPERATING MANUAL BEFORE ATTEMPTING TO USE THE INSTRUMENT

To comply with the standards IEC 60601-1 for safety and IEC 60601-1-2 for EMC the audiometer is designed to be used only with the medically-approved mains adapter supplied, which is specified as part of the equipment. **Do not use any other type of mains adapter with this instrument. Refer to Section 12 for the stock number of the adapter.**

The audiometer is for indoor use only and should be used only as described in this manual.

The transducers supplied with the audiometer are specifically calibrated with it; if these transducers are changed calibration will be required.



Do not immerse the unit in any fluids. See Section 8 of this manual for the proper cleaning procedure for the instrument and its accessories.

Do not use the instrument in an oxygen-rich environment or in the presence of a flammable anaesthetic mixture or other flammable agents.

Do not drop or otherwise impact this instrument. If the instrument is dropped or damaged, return it to the manufacturer for repair and/or calibration. Do not use the instrument if any damage is suspected.

The instrument must be stored and used within the specified temperature, pressure and humidity ranges (see Sections 7 and 9).

Do not attempt to open, modify or service the instrument. Return the instrument to the manufacturer or distributor for all repair and servicing requirements. Opening the instrument will void the warranty.

2.2 Electromagnetic compatibility (EMC) considerations

Medical electrical equipment needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in Appendix 1. This provides guidance on the electromagnetic environment in which to operate the instrument.

Portable and mobile radio-frequency (RF) communications equipment can affect medical electrical equipment. The instrument should not be used adjacent to or stacked with other equipment; if this is necessary the instrument should be observed to verify normal operation.

2.3 Power supply options

The audiometer is designed for continuous operation and may be powered either by a mains adapter (which is supplied, and specified as part of the equipment) or optional internal batteries.



Battery operation

To fit batteries (if configured for this option), remove the battery compartment cover on the base of the audiometer, fit the 4 x 1.5V 'C' batteries supplied (UK only) according to the indications on the battery holder and replace the battery cover.



Batteries should only be changed outside the patient environment. The operator should not touch the battery connectors and the patient simultaneously.

Note: If using batteries, the instrument will automatically switch off approximately 90 seconds after the last key was pressed in order to save battery power. Any test results will be automatically saved.

The display will show "Low Batt" when the battery voltage is low. It is advisable to change the batteries as soon as this happens. Once the voltage of the batteries is too low to operate the instrument, the message "Replace Battery" will appear. Note that local regulations are likely to cover disposal of used batteries.

Mains operation

All other connections must be made **before** connecting the output lead from the adapter into the POWER input socket on the back of the audiometer. Switch on the mains supply - the indicator on the adapter and the POWER indicator on the audiometer will both illuminate green, showing that the instrument is ready for use.

The output from mains adapter is fitted with electronic circuit protection. In case of overload the adapter will shut down and the indicator will be extinguished. When the fault is cleared the adapter will operate as normal.

The input to the mains adapter is protected with a non-replaceable fuse. If this fails the adapter will not operate.

The mains adapter is the mains disconnect device and therefore the audiometer should be positioned such that easy access to the mains adapter is possible.



If a replacement mains adapter is required, please contact Amplivox or your Amplivox distributor.

2.4 Audiometer connections

All the relevant accessory terminals and connections are labelled to ensure correct identification and connection as follows:-

| Socket Label | Socket Type | Colour Code | Connected Part | Notes |
|---------------|--------------------------|-------------|---|-----------|
| RIGHT LEFT | 6.3mm jack 6.3mm jack | Red Blue | Air conduction headset * | |
| PRINTER | RJ12 socket (6-way) | | Printer * | See 2.5 |
| USB | USB Connector Type B | | Computer (via USB port) | See 2.6 |
| N/A | 6 pin mini DIN | | Reserved port; Amplivox diagnostic use only | See below |
| POWER | 2.5mm power jack | | Mains AC/DC Adapter * | |
| RESPONSE | 6.3mm jack | Black | Patient Response Switch * | |

The relevant part numbers are indicated in Section 12

Note regarding the 6-pin mini DIN connector:

This is a restricted socket for Amplivox use only. No user access is permitted.



For connected parts marked * only connect the accessories supplied with the instrument or supplied by Amplivox or an Amplivox distributor. These parts have been tested for use with the Model 116 Screening Audiometer for compliance with the standards IEC 60601-1 and IEC 60601-1-2. The use of accessories other than those specified may compromise compliance with these standards. For other sockets refer to Appendix 2.



2.5 Data transfer to a printer



Please refer to Appendix 2 for important information regarding the connection of non-medical electrical equipment to medical electrical equipment

The audiometer can be upgraded with an option to allow connection to one of three designated portable thermal printers for printing audiometric test results (see Section 3.6). You must use the designated cable for each printer, which is supplied with this option.

Upon receipt of the printer it must be initially charged for a minimum of 15 hours prior to use.

2.6 Data transfer to a computer




Please refer to Appendix 2 for important information regarding the connection of non-medical electrical equipment to medical electrical equipment


The audiometer is supplied with software to allow connection to a computer for the transfer of test results (see Section 3.7). You must use the designated USB cable which is available from Amplivox (see Section 12).



3 Using the Audiometer

3.1 Switching the audiometer on and off

Press and briefly hold the switch marked  (located on the back panel). No warm-up time is required. The display will briefly show the model and the type of headphone currently in use. The display will then be similar to that shown in Section 3.3.

To switch off, press the switch marked  again, or press and hold the MENU key followed by the YES (RIGHT) key and then release both.

3.2 Testing the patient response switch

Press the patient response switch and the light labelled RESPONSE (above and to the right of the display) will illuminate green.

3.3 Audiometer display

On start-up the display will show the following default setting:-

| SIGNAL dBHL | FREQUENCY Hz |
|-------------|--------------|
| 30dB | 1kHz |
| < > | < > |

This indicates that when the PRESENT key is pressed, a tone will be presented at 30dBHL at a frequency of 1kHz (1000Hz) to the designated ear. On start up the audiometer defaults to the left ear.

3.4 Audiometer controls

3.4.1 Multifunction Keys

Several keys on the audiometer have different functions depending on the actual mode of operation. These are MENU (OFF), LEFT (NO), RIGHT (YES) and FREQUENCY \leftrightarrow (MENU SELECT). The use of these keys is described below.



3.4.2 MENU

Press and hold MENU to access the following options. Use the MENU SELECT keys to step through the available options and then the NO, YES or SIGNAL ↓↑ keys to select an action or modify a setting. Release of the MENU key then initiates the action or saves the modified setting and returns to the default display.

| <u>Menu Option</u> | <u>Description</u> |
|--------------------------|---|
| Switch off?: | As described in Section 3.1 |
| Clear test?: | Press YES and release MENU to clear the Threshold Retention Function results from the previous test (see also Section 3.5) |
| Contrast: | Adjust contrast using the SIGNAL ↓↑ keys |
| Store on 2 of 3?: | If activated a hearing threshold will be stored automatically when the patient makes a response to 2 out of 3 tone presents (see Section 3.5.2) |
| Default level: | Adjust the default tone presentation level using the SIGNAL ↓↑ keys |
| Select printer: | Use the SIGNAL ↓↑ keys to select either the Able AP1300 or the Sanibel MPT-II printer |

3.4.3 Description of Function of Other Keys

| | |
|--------------|---|
| PRINT | Press to print the threshold levels displayed (see Section 3.6) |
| LEFT | Press once to select the left ear (the indicator above the key illuminates green); if the left ear is already selected press again to store the displayed signal value as a threshold (see Section 3.5.1) |



| | |
|-----------------------|--|
| RIGHT | Press once to select the right ear (the indicator above the key illuminates green); if the right ear is already selected press again to store the displayed signal value as a threshold (see Section 3.5.1) |
| SIGNAL | Press the ↓↑ keys to decrease or increase the level of the tone presented in 5dB steps; to scroll through the range keep the key pressed |
| PULSE/CONSTANT | Press this key once to enable the pulse tone present function when the PRESENT key is operated (“Pulse” is displayed); press the key again to enable constant tone presentation (the indicator above the key illuminates green) and use the PRESENT key to interrupt the tone; press the key once again to return to default |
| +20dB | This enables tone levels to be presented with up to 20dB higher output; press the key and then use SIGNAL ↑ to access the extra 20dB in 5dB steps; an indicator above the key illuminates green to show that the function is active Note: this function will not operate if constant tone presentation has been enabled |
| WARBLE | When selected, the instrument presents a warble tone a green indicator above the key is used to show that the function is active; this feature may be used in conjunction with the PULSE/CONSTANT key |
| FREQUENCY | Press the ⇐ key to select a lower frequency and the ⇒ key to select a higher frequency |
| TALK OVER | Hold this key to interrupt the test and route the operator’s voice from the front panel |



microphone to the headset; the level is adjusted with the SIGNAL ↓ ↑ keys; when the key is released the output levels return to their default values

PRESENT

Press to present the displayed test signal to the patient. The “PRESENT” indicator above the display will be illuminated green during tone presentation.

3.5 Threshold Retention Function

This function records the thresholds for both ears at each frequency tested. Thresholds may be recorded manually or automatically.

The operator can then review the results at the end of the test and record them on an audiogram card, print them with the optional printer (see Section 3.6), or transfer the results to a computer (see Section 3.7).

3.5.1 Recording thresholds manually

Once a threshold has been determined press the “selected” ear key once again. The threshold will be recorded and displayed as shown in the illustration in Section 3.5.3. **Note: this function will not operate if the “Store on 2 of 3” option has been enabled (see Section 3.5.2).**

3.5.2 Recording thresholds automatically

If the “Store on 2 of 3” option has been enabled (see Section 3.4.2) then a threshold will be recorded automatically by the audiometer if the patient makes a response to two out of three manual tone presentations at the same level and frequency. Thresholds determined using the “Store on 2 of 3” option are displayed within square brackets.



3.5.3 Reviewing retained thresholds

To review the retained thresholds, select the required frequency using the FREQUENCY $\leftarrow \rightarrow$ keys. The recorded values for the left and right ears are shown on the lower line of the display as illustrated below.

| SIGNAL dBHL | FREQUENCY Hz |
|-------------|--------------|
| 30dB | 4kHz |
| 20 | 10 |

THRESHOLDS

This display shows thresholds at 4kHz
Left ear 20dBHL
Right ear 10dBHL

To clear the Threshold Retention memory, use the Clear Test menu option described in Section 3.4.2.

3.6 Printing audiograms

Two designated thermal printers (the Able AP1300 or the Sanibel MPT-II) are available as options for use with the Model 116 audiometer. The correct printer must be selected (use the MENU options described in Section 3.4.2 to make this selection).

- Connect the PRINTER socket on the audiometer (6-way RJ12) to the printer with the supplied cable (refer to Section 2.5 of this operating manual for printer set-up). **Note that the printer cables for the Able printer (A108) and Sanibel (A102) are not compatible.**
- Ensure that the printer is fully charged, switched on, loaded with paper and ready to print.
- Press the PRINT key and on the prompt "Is printer ready?" press the YES key. The audiogram will then print. To cancel press NO.

3.7 Data transfer to Audibase or ampliSuite

Test results stored within the audiometer may be transferred to the Amplivox Audibase database which is available as an option and must be installed on to a computer (see Section 12 for the part number). Alternatively, Amplivox ampliSuite allows data to be transferred to a computer and subsequently viewed, annotated & printed. This software is supplied on a CD which includes this operating manual.



Refer to the installation & operating instructions provided with Audibase or ampliSuite for further details.

4 Suggested Sequence of Operation and Test Procedure

The following applies to air conduction measurements. Refer also to ISO 8253 for guidance.


4.1 Audiometry preparation and ambient conditions

Refer to the appropriate audiometric standards and other relevant publications for guidance on audiometric testing.

Audiometric testing should always be performed in quiet conditions (e.g. a quiet room or an acoustic booth). The optional Audiocups can provide an additional level of isolation from ambient noise. For further explanation on permissible ambient noise levels, please refer to the standard ISO6189.

4.2 Test system arrangement

The schematic diagram below shows a typical example of the use of audiometric test equipment. The audiometer is located on the desk of a seated operator as shown.

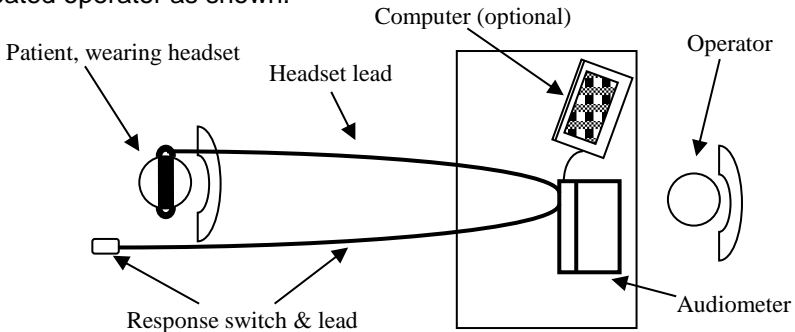

The patient is seated in front of the desk facing away from the operator. The patient wears a headset (see Section 4.3) and responds to test stimuli by use of a hand-held switch which is also connected to the instrument.



4.3 Headset

The headset must be fitted by a qualified person to ensure a proper seal and a comfortable fit. The leads from the headset are connected to the instrument and the headset is then fitted to the patient.

4.4 Patient instructions

The patient should be given the following instructions using the TALKOVER function:

- “As soon as you hear the tone, press the response switch. When you no longer hear the tone release the response switch”

4.5 Pre-test

- (1) Switch the audiometer on
- (2) Perform a listening check
- (3) Decide whether to use the manual or automatic Threshold Retention Function and/or an audiogram card to record the thresholds
- (4) If the automatic Threshold Retention Function is required ensure that the Store on 2 of 3 option is enabled (see Section 3.5.2) and that a patient response switch is in use
- (5) Prepare the test environment & patient (see Sections 4.1 to 4.4)
- (6) If the patient response switch is not being used give instructions to the patient to acknowledge any tone presented by raising or lowering the finger
- (7) Select the better hearing ear (according to the patient) by pressing either the LEFT or RIGHT key and start the familiarisation session.

4.6 Familiarisation

- (1) Present the tone 30dB at 1kHz for between 1 and 2 seconds. If there is no response at 30dB, increase the attenuation level in 10dB steps until the patient responds
- (2) When the patient responds, wait for 1 to 2 seconds and present the tone again at the same level; however, if the patient does respond at 30dB, reduce the signal level in 10dB steps, repeating the presentation until there is no response, then increase the signal level



in 5dB steps until the patient responds; wait 1 to 2 seconds and present the tone again at the same level

- (3) If the responses are consistent with the pattern of tone presentation proceed to Section 4.7 and start measuring the patient's hearing thresholds; if not, repeat the familiarisation process

4.7 Test

- (1) Present the first test tone at 30dB at 1kHz
- (2) If the patient responds, reduce the signal level in 10dB steps repeating the presentation until there is no response; then increase the signal level in 5dB steps until the patient responds
- (3) If the patient does not respond, increase the signal level in 5dB steps until there is a response and then continue with step 4.
- (4) Repeat the test by reducing the signal level in 10dB steps until the patient no longer responds. Then increase the signal level in 5dB steps until they do respond and note this level.
- (5) Repeat step 4 until the patient responds three out of a maximum of five times at the same signal level. This indicates the patient's hearing threshold level for that frequency. Either mark the threshold on an audiogram card or press the appropriate ear key once to activate the Threshold Retention Function and save the threshold level on screen.
- (6) Proceed to the next test frequency. It is common practice to test the frequencies in the following order: 1k, 2k, 3k, 4k, 6k, 8k and 500 Hz.
- (7) Repeat steps 1 to 6 for the other ear.

4.8 Post-test

- (1) Use the Threshold Retention Function to review the results (See 3.5)
- (2) If required do one or more of the following:
 - Record the results on an audiogram card, or
 - Print the results (Section 3.6), or
 - Transfer the results to a computer (Section 3.7)

Refer to Section 3.4.2 to clear the thresholds at the end of a test and, if required, switch off the audiometer.



5 Specification

5.1 Output data

| | |
|-------------------------|--|
| Outputs: | Left and Right earphone |
| Frequency range: | 125Hz to 8kHz |
| Frequency accuracy: | <1% |
| Distortion: | <2% |
| Output level range: | -10dBHL min; see Section 5.2 for maximum |
| Output level accuracy: | Within 3dB |
| Output level step size: | 5dB |
| Output transducer: | DD45 earphones (supplied) |
| Tone present: | Single, warble or pulsed |
| Communication: | Integral talk over facility |
| USB interface: | Transfer of test results to a computer |

5.2 Maximum hearing levels provided at each frequency

| Frequency, Hz | Air conduction, dBHL (DD45) | Frequency, Hz | Air conduction, dBHL (DD45) |
|---------------|-----------------------------|---------------|-----------------------------|
| 125 | 70 | 2000 | 100 |
| 250 | 80 | 3000 | 100 |
| 500 | 90 | 4000 | 100 |
| 750 | 100 | 6000 | 100 |
| 1000 | 100 | 8000 | 80 |
| 1500 | 100 | | |

5.3 Physical Data

| | |
|---------------------------|--|
| Display: | 2 lines of 24 characters |
| Battery power (optional): | 4x1.5V "C" cells (alkaline recommended) |
| Mains Power: | 100-240Vac; 50-60Hz; 0.5A |
| Input Rating: | 5Vdc; 1.2 A |
| Dimensions: | 270mm long x 165mm deep x 60mm high |
| Weight (no batteries): | 0.75kg (approx) |
| Safety: | IEC 60601-1 (plus UL, CSA & EN deviations) |
| EMC: | IEC 60601-1-2 |
| CE mark: | To the EU Medical Device Directive |



5.4 Equipment classification

| | |
|---|---|
| Type of protection against electric shock | Powered via SELV ClassII mains adapter |
| Degree of protection against electric shock | Type B applied part |
| Degree of protection against ingress of water | Not protected |
| Mode of operation | Continuous operation |
| Equipment mobility | Portable |

The Model 116 Audiometer is classified as a Class IIa device under Annex IX of the EU Medical Devices Directive. It is intended for use as a screening audiometer instrument.



6 Symbols

The following symbols appear on the audiometer or mains adapter:



Definition: Identifies the control by means of which the instrument is switched on from (or returned to) a standby condition.



Definition: Refer to instruction manual (mandatory).




Definition: Type B applied part – an applied part providing protection against electric shock, particularly regarding allowable patient leakage current and patient auxiliary current.

The applied parts are the left & right earphones, patient response switch and the associated cables.



Definition: The output from the mains AC adapter is Direct Current.

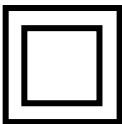

Definition: Class II equipment – equipment in which protection against electric shock does not rely on basic insulation only, but in which additional safety precautions such as double insulation or reinforced insulation are provided, there being no provision for protective earthing or reliance upon installation conditions.



7 Technical Information

Audiometer

Audiometer type: Type 4 (IEC 60645-1:2001)
Type 4 (ANSI S3.6:2004)

Battery function

Battery voltage range: 4.0 to 6.0V
Low battery warning: Approx 4.4V
Expected battery life: 6 to 8 hours use from alkaline batteries

Frequency Modulation

Carrier frequencies: 125Hz to 8kHz as per pure tones
Modulation waveform: Sinusoidal
Rising and falling symmetry: Symmetrical on linear frequency scale
Modulating frequency: 15.625Hz
Frequency deviation: +/-10%

Transducers

Types and reference levels: DD45: ISO 389-1, Table 2
Static headband force: Headphones: 4.5N
Sound attenuation characteristics: ISO8253-1, Table 3

Earphone Sound Attenuation Characteristics

| Frequency, Hz | 125 | 250 | 500 | 1000 | 2000 | 4000 | 8000 |
|-----------------|-----|-----|-----|------|------|------|------|
| Attenuation, dB | 2 | 5 | 7 | 15 | 25 | 31 | 23 |

Environmental

Operating temperature: +15°C to +35°C
Operating humidity: 30% to 90% (non-condensing)
Atmospheric pressure: 700 hPa to 1060 hPa

Input / Output

Power input: 2.5mm barrel-type socket.
Patient response input: 6.3mm Jack socket
Left & Right outputs: 6.3mm Jack socket
USB: Type B socket
Printer: RJ12 socket (6-way)
Maximum voltage at any output: 12V peak



8 Routine Maintenance

8.1 Audiometer maintenance

The Model 116 audiometer is a precision instrument. Handle it carefully in order to ensure its continued accuracy and service. When cleaning the instrument, first disconnect it from the mains supply. Use a soft cloth and mild detergent to clean the instrument panel when required. Refer to ISO 8253-1 for additional guidance.

8.2 Transducer maintenance

Before use check the transducer cables and connectors for signs of wear and/or damage. If you find any, please replace the item immediately by contacting Amplivox or your Amplivox distributor, requesting the relevant part number (see Section 12).

Handle the audiometric headset and other accessories with care. For parts that are in direct contact with the patient it is recommended that replacement parts are used or the parts are subjected to a standard disinfecting procedure between patients.

This includes physically cleaning and use of a recognised disinfectant. The specific manufacturer's instructions should be followed for use of this disinfecting agent to provide an appropriate level of cleanliness.

Clean the ear cushions (including those on the Audiocups, if used) with a recognised disinfectant, e.g. a "Mediswab".



During the cleaning process do not allow moisture to enter the earphone.



8.3 Mains adapter maintenance

Before use check the mains AC adapter for signs of wear and/or damage. If you find any replace the adapter immediately by contacting Amplivox or your Amplivox distributor. Refer to Section 12 for approved part numbers.



DO NOT USE ANY OTHER TYPE OF MAINS ADAPTER WITH THIS INSTRUMENT. See Section 2.3.

8.4 Batteries

Batteries (if fitted) should be removed if the instrument is not to be used over an extended period of time.

9 Instrument Storage and Transportation

This instrument can be stored or transported within the following environmental parameters:

| | |
|-----------------------|-----------------------------|
| Temperature: | -20°C to +70°C |
| Humidity: | 10% to 90% (non-condensing) |
| Atmospheric Pressure: | 500 hPa to 1060 hPa |



10 Calibration and Repair of the Instrument

Amplivox recommends that this audiometer should be calibrated on an annual basis. Please contact Amplivox or the designated distributor for details of calibration services. Refer to ISO 8253-1 for additional guidance.



The instrument should be returned to the manufacturer for service & repair. There are no user-serviceable parts within it.

When packing the instrument for shipping, please use the original shipping carton and packing materials. Please also ensure that the headset leads are not wrapped around the headband of the headset.

11 Guarantee

All Amplivox instruments are guaranteed against faulty materials and manufacture. The instrument will be repaired free of charge for a period of two years from the date of despatch if returned, carriage paid, to the Amplivox service department. Return carriage is free of charge for customers in the UK and chargeable for overseas customers.

Important Note:

The following exceptions apply:

Earphones may go out of calibration due to rough handling or impact (dropping). The life of the leads is also dependent upon conditions of use. These parts are only guaranteed against faulty materials or manufacture.



12 Ordering Consumables and Accessories

To order consumables, additional accessories and to replace detachable parts that have been damaged, please contact Amplivox for current prices and delivery charges. The items available are listed below:

| Stock No. | | Description |
|-----------|---------|--|
| A022 | 8010855 | Audiocups (noise reducing earphone enclosures) |
| AC1042 | 8010835 | Audiocup ear cushion |
| AC1047 | 8507381 | Audiocup headband |
| AC1048 | 8010834 | Audiocup headband cover |
| | | |
| A023 | 8010840 | Headband (standard headphone) |
| A026 | 8010857 | Earphone cushion |
| A032 | 8010876 | Earphones DD45 * |
| A030 | 8010822 | Headset lead |
| | | |
| B128 | 8004674 | Carrying case |
| A091-7 | 8512734 | Approved mains adapter |
| A085 | 8103838 | Patient response switch |
| A051 | 8013007 | Audiogram cards (pack of 50) |
| | | |
| PT02 | 8029305 | Printer Sanibel MPT-II |
| A102 | 8004419 | Printer cable for audiometer to Sanibel MPT-II |
| C0104 | 8029305 | Thermal Printer paper for Sanibel MPT-II |
| | | |
| F07 | 8507230 | USB Cable, 2.0m |
| AUD06 | 8511500 | Amplivox Audibase 5.5 (including USB cable) |



Accessories marked * require calibration with the specific audiometer to be used. Do not attempt to use these accessories until the audiometer has been calibrated to match their characteristics.

Shipping documentation will reference the stock number quoted above, and images of the parts alongside the relevant stock number are available



on the Amplivox website (www.amplivox.ltd.uk). The required fitting instructions are supplied with each part.

13 Disposal Information



Amplivox Limited is fully compliant with the WEEE (Waste Electrical and Electronic Equipment) Regulations. Our PRN (Producer Registration Number) is WEE/GA0116XU and we are registered with the approved WEEE Compliance Scheme, B2B Compliance, approval number WEE/MP3338PT/SCH.

The main purpose of the WEEE Regulations is to encourage the segregation of waste electrical items from the general waste stream and into reuse, recovery and recycling routes.

For any waste electrical units purchased from Amplivox that either:

- bear the crossed out wheeled bin symbol with black bar underneath
- or, have been replaced with new Amplivox products on a like-for-like basis

please contact our WEEE Compliance Scheme using the details below. B2B Compliance will be able to provide further information on how to recycle your waste electrical units and answer any queries you may have.

B2B Compliance

Tel: +44 (0) 1691 676 124 (Option 2)

Email: operations@b2bcompliance.org.uk



Appendix 1 - EMC Guidance & Manufacturer's Declaration

| Guidance and manufacturer's declaration – electromagnetic emissions | | |
|---|-------------------|---|
| The Model 116 Audiometer is intended for use in the electromagnetic environment specified below. The customer or user of Model 116 Audiometer should assure that it is used in such an environment. | | |
| Emissions test | Compliance | Electromagnetic environment – guidance |
| RF emissions CISPR 11 | Group 1 | The Model 116 Audiometer uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause interference in nearby electronic equipment. |
| RF emissions CISPR 11 | Class B | The Model 116 Audiometer is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes. |
| Harmonic emissions IEC 61000-3-2 | Class A | |
| Voltage fluctuations/flicker emissions IEC 61000-3-3 | Complies | |



Guidance and manufacturer's declaration – electromagnetic immunity (1)

The Model 116 Audiometer is intended for use in the electromagnetic environment specified below. The customer or user of the Model 116 Audiometer should assure that it is used in such an environment.

| Immunity test | IEC 60601 test level | Compliance level | Electromagnetic environment – guidance |
|--|--|--|--|
| Electrostatic Discharge (ESD) IEC 61000-4-2 | ±8 kV contact ±15 kV air | ±8 kV contact ±15 kV air | Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30% |
| Electrical fast transient/burst IEC 61000-4-4 | ±2 kV for power supply lines ±1 kV for input/output lines | ±2 kV for power supply lines ±1 kV for input/output lines | Mains power quality should be that of a typical commercial or hospital environment |
| Surge IEC 61000-4-5 | ±1 kV differential mode ±2 kV common mode | ±1 kV differential mode ±2 kV common mode | Mains power quality should be that of a typical commercial or hospital environment |




| Immunity test | IEC 60601 test level | Compliance level | Electromagnetic environment – guidance |
|---|--|--|---|
| Voltage dips, short interruptions and voltage variations on power supply input lines (100V/60Hz & 240V/50Hz) IEC 61000-4-2 | 0% UT (100% dip in UT) for 0.5 cycle 0% UT (100% dip in UT) for 1 cycle 40% U_T (60% dip in U_T) for 5 cycles 70% UT (30% dip in UT) for 500ms 0% UT (100% dip in UT) for 5 sec | 0% UT (100% dip in UT) for 0.5 cycle 0% UT (100% dip in UT) for 1 cycle 40% U_T (60% dip in U_T) for 5 cycles 70% UT (30% dip in UT) for 500ms 0% UT (100% dip in UT) for 5 sec | Mains power quality should be that of a typical commercial or hospital environment. If the user of the Model 116 Audiometer requires continued operation during power mains interruptions, it is recommended that the Model 116 Audiometer be powered from an uninterruptible power supply or a battery |
| Power frequency (50/60 Hz) magnetic field IEC 61000-4-8 | 30 A/m | 30 A/m | Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment. |
| NOTE U_T is the a.c. mains voltage prior to the application of the test level | | | |



| Guidance and manufacturer's declaration – electromagnetic immunity (2) | | | |
|---|-------------------------------|-------------------------------|---|
| The Model 116 Audiometer is intended for use in the electromagnetic environment specified below. The customer or user of the Model 116 Audiometer should assure that it is used in such an environment. | | | |
| Immunity test | IEC 60601 test level | Compliance level | Electromagnetic environment – guidance |
| Conducted RF IEC 61000-4-6 | 10 Vrms 150kHz to 80MHz | 10 Vrms 150kHz to 80MHz | Portable and mobile RF communications equipment should be used no closer to any part of the Model 116 Audiometer, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 1.2\sqrt{P}$ |
| Radiated RF IEC 61000-4-3 | 10 V/m 80MHz to 2.7GHz | 10 V/m 80MHz to 2.7GHz | $d = 1.2\sqrt{P}$ 80MHz to 800MHz $d = 2.3\sqrt{P}$ 800MHz to 2.5GHz where P is the maximum output power rating of the transmitter in Watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). |



Guidance and manufacturer's declaration – electromagnetic immunity (2)

| | | | |
|--|--|--|--|
| | | | <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^a should be less than the compliance level in each frequency range. ^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p>  |
|--|--|--|--|

NOTE 1 At 80MHz and 800MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

- a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Model 116 Audiometer is used exceeds the applicable RF compliance level above, the Model 116 Audiometer should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Model 116 Audiometer.
- b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.



| Guidance and manufacturer's declaration – electromagnetic immunity (2) | | | |
|--|--|--------------------------|---------------------------|
| Recommended separation distances between portable and mobile RF communications equipment and the Model 116 Audiometer | | | |
| The Model 116 Audiometer is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Model 116 Audiometer can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Model 116 Audiometer as recommended below, according to the maximum output power of the communications equipment. | | | |
| Rated maximum output power of transmitter | Separation distance according to frequency of transmitter | | |
| | m | | |
| W | 150 kHz to 80 MHz | 80 MHz to 800 MHz | 800 MHz to 2.5 GHz |
| | $d = 1.2\sqrt{P}$ | $d = 1.2\sqrt{P}$ | $d = 2.3\sqrt{P}$ |
| 0.01 | 0.12 | 0.12 | 0.23 |
| 0.1 | 0.38 | 0.38 | 0.73 |
| 1 | 1.2 | 1.2 | 2.3 |
| 10 | 3.8 | 3.8 | 7.3 |
| 100 | 12 | 12 | 23 |



Guidance and manufacturer's declaration – electromagnetic immunity (2)

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in Watts (W) according to the transmitter manufacturer.

NOTE 1 At 80MHz and 800MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

NOTE 3 WARNING: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30cm (12 inches) to any part of the Model 116 audiometer including cables specified by the manufacturer. Otherwise degradation of the performance of this equipment could result.



Appendix 2 - Use with Non-medical Electrical Equipment

Any person who connects external equipment to signal input, signal output or other connectors has created a medical electrical system and is therefore responsible for the system complying with the requirements of clause 16 of IEC 60601-1:2005 (*General requirements for basic safety and essential performance*).

If connections are made to standard equipment such as printers and computers, special precautions must be taken in order to maintain medical safety. The following notes are provided for guidance in making such connections to ensure that the general requirements of clause 16 of IEC 60601-1:2005 are met.

The following signal inputs and outputs on the Model 116 audiometer are electrically isolated to the requirements of IEC 60601-1 in order to reduce any potential hazard associated with the use of mains-powered equipment connected to these inputs and outputs:

| Socket Label | Socket Type | Typical Connection |
|--------------|----------------------|--------------------|
| PRINTER | RJ12 socket (6-way) | Printer |
| USB | USB Connector Type B | Computer |

External equipment intended for connection to signal input, signal output or other connectors, shall comply with the relevant IEC or international standards (e.g. IEC 60950, CISPR 22 & CISPR 24 for IT equipment, and the IEC 60601 series for medical electrical equipment).

Equipment not complying with IEC 60601 shall be kept outside the patient environment, as defined in IEC 60601-1:2005 (at least 1.5m from the patient).

The operator must not touch the connected equipment and the patient at the same time as this would result in an unacceptable hazard.



Refer to Diagrams 1 to 3 below for typical configurations of connected peripheral equipment. Refer to Amplivox Limited at the address given on the front of this user manual if advice is required regarding the use of peripheral equipment.

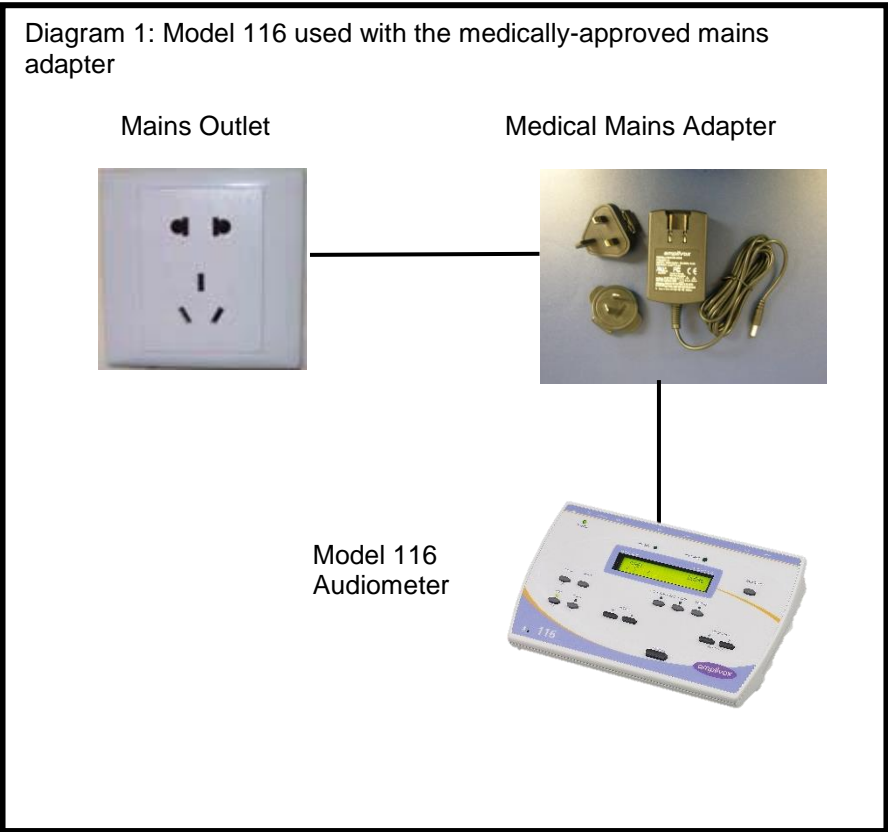


Diagram 2: Model 116 used with the medically-approved mains adapter and printer


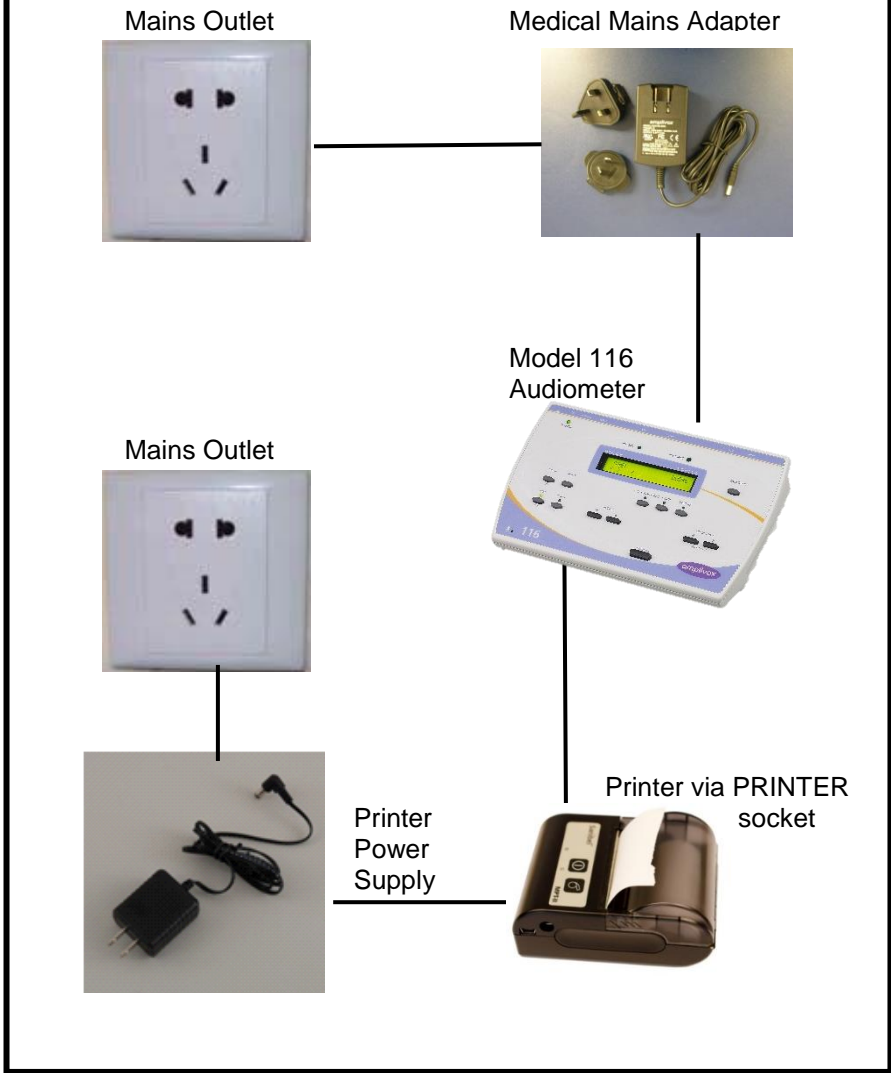


Diagram 3: Model 116 used with the medically-approved mains adapter and PC



Copyright © 2018 Amplivox Ltd.

All rights reserved. No part of this publication may be reproduced or transmitted in any form or by any means without the prior written permission of Amplivox Ltd.

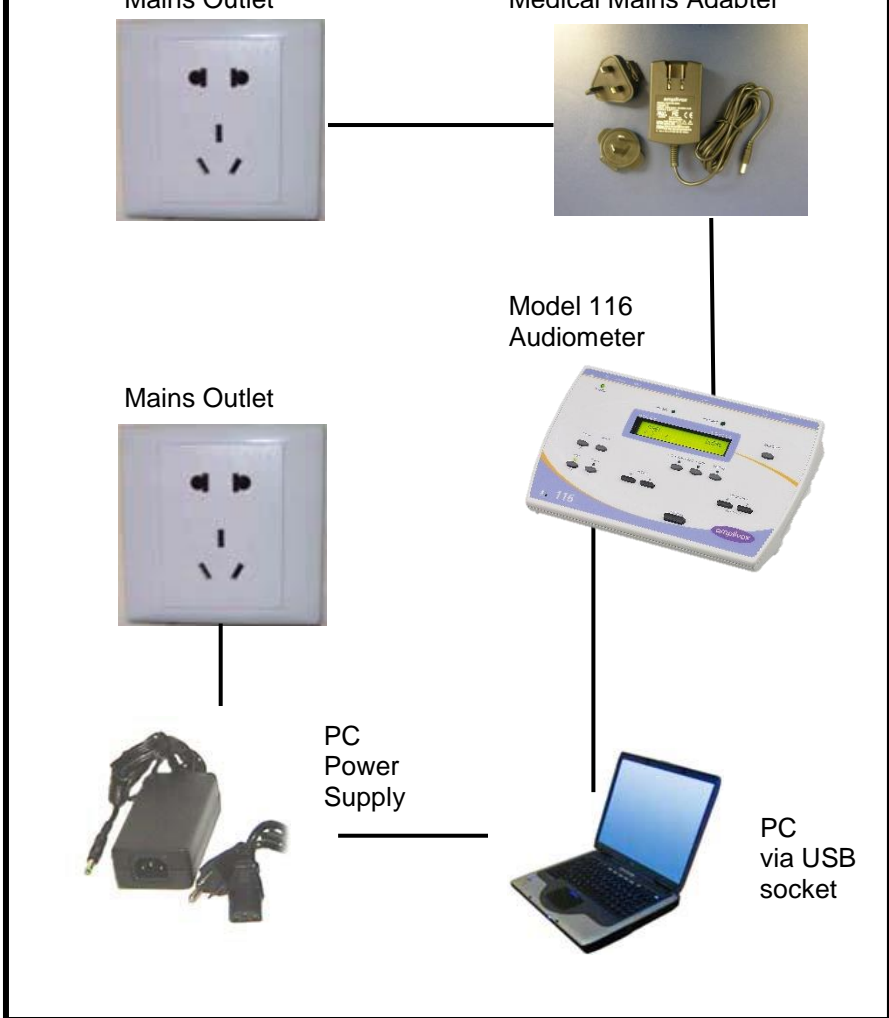