

Instructions for use



Charger 1.0 miniRITE R

Oticon Opn S™
Oticon Opn Play™
Oticon Ruby

Introduction to this booklet

This booklet guides you on how to use and maintain your new charger. Please read the booklet carefully, including the **Warnings** section. This will help you get the most out of the charger. If you have additional questions about the charger, please contact your hearing care professional.

| **About** | [Start-up](#) | [Handling](#) | [Warnings](#) | [More info](#) |

For your convenience, this booklet contains a navigation bar to help you easily navigate through the different sections.

Intended use

Intended use	The charger is intended to recharge hearing aids that contain rechargeable batteries. The product is intended to be used together with the given hearing solution.
Indications for use	There are no indications for use (diagnoses) for the Charger itself.
Intended user	Lay persons including children older than 36 months.
User environment	Home setting.
Contraindications	No contraindications.
Clinical benefits	See clinical benefits of the hearing aid.

IMPORTANT NOTICE

Please refer to your hearing aid's instructions for use for details on the functionality of your specific hearing aid.

Table of content

About

Charger	6
Temperature conditions	7

Start up

Prepare charger	8
Charge the hearing aids	9
Charging time	10

Handling

Caring for your charger	12
Power the charger using other sources	13

Warnings


Warnings	14
----------	----

More info

Troubleshooting	16
Charger: Conditions of use	20
Warranty	21

Charger

What it is and does

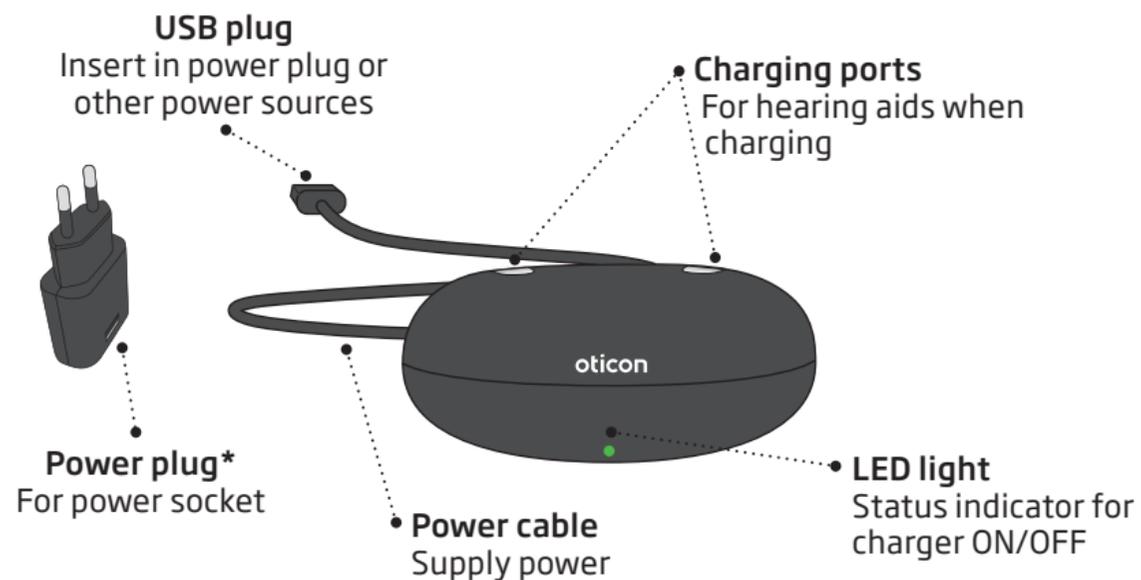

* Power plug will vary from country to country


Temperature conditions

Temperature	Charging condition
Below +5°C	No charging due to low temperature
+10°C to +35°C	Best charging condition. Charging will take up to 3 hours.
+5°C to +10°C and +35°C to +40°C	Prolonged charging. Charging will take up to 4 hours.
Above +40°C	No charging due to high temperature

Charging will automatically resume when the temperature is within the operating range of +5°C to +40°C.

Prepare charger


Connect to power source

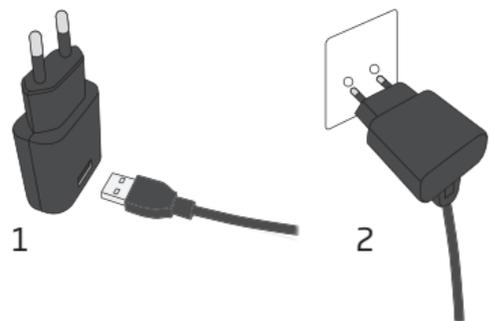

1. Insert the USB plug into the power plug.
2. Insert the power plug into a power socket.

The charger will turn on automatically.

Charger is on

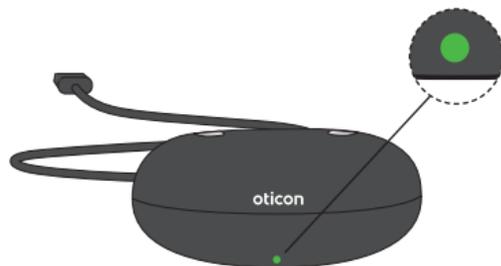

When the charger is connected to power, the green LED turns on.

Charge the hearing aids


Place hearing aids in charger



Place the hearing aids in the charging ports as shown in the picture. Make sure that the hearing aids match the right and left indicators on the charger.

The right hearing aid should be placed where the red indicator is. The left hearing aid should be placed where the blue indicator is.

Charging

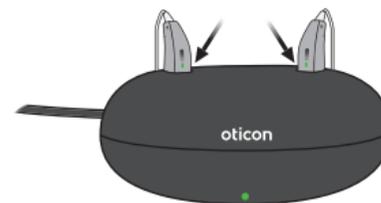

The hearing aids LEDs will turn on when the hearing aids are correctly placed in the charger.

Hearing aid LED charging status:
Red /Orange* = Charging
Green light = Fully charged

*depending on the hearing aid family, see hearing aid instruction for use.

Charging time

Ensure you fully charge your hearing aid before first time use.

Charge your hearing aid every night. That ensures you start your day with a fully charged hearing aid.

Charging time may vary depending on the remaining capacity of the battery. If your hearing aid is completely drained the normal charging time is:

3 hours	1 hour	0.5 hour
Fully charged	50% charged	25 % charged

The charging process will automatically stop when the battery is fully charged. Leave the hearing aid in the charger when you do not use it. Remember to keep the charger connected to a power source.

Note charging time may vary between left and right hearing aid.

IMPORTANT NOTICE

Always leave the charger connected to a power source when the hearing aid is seated in the charging port. Switching off the charger will make the hearing aid turn on and start using battery.

Caring for your charger

Always make sure that the charger is dry and clean before use. Keep the charging ports clean for dust using a dry cloth. Do not use water or cleaning fluids to clean the charging ports.

Store your charger in a dry location, e.g. in the bedroom, rather than in the bathroom or kitchen.

Keep the charger in a temperature environment of +5°C to +40°C.

Power the charger using other sources

It is possible to power the charger from other sources with a USB port. Ensure that the power source is USB 2.0 compliant, minimum 500mA output.

Examples of sources:

- Power bank
- PC
- Car

When using a power bank ensure that it does not have a power-saving or auto-shutdown mode. As this can result in your hearing aids not being fully charged or not charged at all.

IMPORTANT NOTICE

When recharging the hearing aid use the supplied original power supply 5.0V/1 A, CE and safety marked, for the charger.



Warnings

For your personal safety and to ensure correct usage, you should familiarise yourself fully with the following general warnings before using your charger. Consult your hearing care professional if you experience unexpected operations or serious incidents with your charger during use or because of its use. Your hearing care professional will support you with issue handling and, if relevant, reporting to the manufacturer and/or the national authorities.

Usage of charger

The charger should be used only as directed and adjusted by your hearing care professional.

Risk of burns

Do not place anything but the hearing aids in the charger ports.

Choking hazards

Charger, hearing aids and their parts should be kept out of reach of children and anyone who might swallow these items or otherwise cause injury to themselves.

If a small part is swallowed, see a doctor or hospital immediately.

Heat

The charger must never be exposed to extreme heat, e.g. left inside a parked car in the sun.

The charger must not be dried in microwave ovens or other ovens.

Connection to external equipment

The safety of the use charger when connected to external equipment with an auxiliary input cable and/or with USB cable and/or directly, is determined by the external signal source. When connected to external equipment plugged into a wall outlet, this equipment must comply with IEC-62368-1 or equivalent safety standards.

Fire hazards

The power source in your charger has insufficient energy to cause fire in normal usage conditions.

Troubleshooting

Symptom	Possible causes	Solutions
The charger LED remains turned off when the charger is connected to power	The charger is not turned on	Verify that the power supply is connected correctly Always control charging starts up to make sure the power source delivers sufficient output for your charger
The charger LED shows long orange blinks	The hearing aid or charger is too warm	Move the charger and hearing aid to a location with a temperature between +5°C to +40°C
The charger LED repeatedly shows 3 short orange blinks	The charger is too cold	Move the charger to a location with a temperature between +5°C to +40°C
	System error	Contact your hearing care professional

Troubleshooting

Symptom	Possible causes	Solutions
Hearing aid LED remains off when the hearing aid is placed in the charger	The charger is not turned on	Verify that the charger power supply is connected correctly
	The hearing aid or charger is either too warm or too cold	Move the charger and hearing aid to a location with a temperature between +5 °C to +40 °C
	Charging is incomplete. The room temperature exceeds +35 °C, which prolongs the charging time. The charger has stopped charging to protect the battery.	Reinsert the hearing aid into the charger. This will complete the charging within approximately 15 minutes
	The hearing aid is not correctly inserted into the charger	Check the charger slots for foreign objects
The hearing aid LED blinks red/orange* when the hearing aid is placed in the charger	System error	Contact your hearing care professional

* depending on the hearing aid family, see hearing aid instruction for use

Charger: Conditions of use

Operational environment conditions	Temperature: +5°C to +40°C Humidity: 5% to 93%, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa
Transportation and storage conditions	Temperature and humidity shall not exceed the below limits for extended periods during transportation and storage: Temperature: -25°C to +70°C Humidity: 5% to 93%, non-condensing

Information about hearing aid: Condition of use, see hearing aid Instruction for use.

Warranty

Certificate

Name of owner: _____

Hearing care professional: _____

Hearing care professional's address: _____

Hearing care professional's phone: _____

Purchase date: _____

Warranty period: _____ Month: _____

Model: _____ Serial no.: _____

International warranty

Your charger is covered by an international limited warranty issued by the manufacturer for a period of 12 months from the date of delivery. This limited warranty covers manufacturing and material defects in the device itself, but not in accessories such as e.g. batteries, cables, power supplies, etc. Problems arising from improper/incorrect handling or care, excessive use, accidents, repairs made by an unauthorised party, exposure to corrosive conditions, damage due to foreign objects entering the device or incorrect adjustments are NOT covered by the limited warranty and may void it. The above warranty does not affect any legal rights that you might have under applicable national legislation governing sale of consumer goods. Your hearing care professional may have issued a warranty

that goes beyond the clauses of this limited warranty. Please consult him/her for further information.

If you need service

Take your charger to your hearing care professional, who may be able to sort out minor problems and adjustments immediately.

The manufacturer declares that this charger is in compliance with the essential requirements and other relevant provisions of Directive EMC (2014/30/EU) and LVD (2014/35/EU).

Declaration of Conformity is available from the manufacturer.



SBO Hearing A/S
Kongebakken 9
DK-2765 Smørum
Denmark
www.sbohearing.com



Waste from electronic equipment must be handled according to local regulations.



Description of symbols and abbreviations used in this booklet	
	Warnings Text marked with a warning symbol must be read before using the device.
	Manufacturer The device is produced by the manufacturer whose name and address are stated next to the symbol. Indicates the medical device manufacturer, as defined in EU Directives 90/385/EEC, 93/42/EEC and 98/79/EC.
	Electronic waste (WEEE) Recycle hearing aids, accessories or batteries according to local regulations. Hearing aid users can also return electronic waste to their hearing care professional for disposal. Electronic equipment covered by Directive 2012/19/EU on waste and electrical equipment (WEEE).
	Regulatory Compliance Mark (RCM) The device complies with electrical safety, EMC and radio spectrum requirements for devices supplied to the Australian or New Zealand market.
	Double insulation Indicates that product power supply is double insulated. This means a very low risk of single fault failures i.e. single fault failures cannot occur.
	Indoor use Indicates that the product is qualified for indoor use.

214213UK / 2019.11.06

Description of additional symbols used on labels	
	Keep dry Indicates a medical device that needs to be protected from moisture.
	Caution symbol Consult instructions for use for warnings and cautions.
	Catalogue number Indicates the manufacturer's catalogue number so that the medical device can be identified.
	Serial number Indicates the manufacturer's serial number so that a specific medical device can be identified.
	Medical Device The device is a medical device.
	Power supply outlet Indicates the rate at which power is sent from the outlet to your electrical device.

214213UK / 2019.11.06



oticon
life-changing technology