

# THE HearingReview

Hearing Loss ▾ Products ▾ Practice Building ▾ Inside Hearing ▾ Resources ▾ Magazine

## How Do Medical Masks Degrade Speech Reception?

Apr 1, 2020 | Assistive Technologies, Counseling & Education, Health & Wellness, Industry News | ★★★★★



Patient Care | May 2020 *Hearing Review*

**Speech blocked by surgical masks becomes a more important issue in the Era of COVID-19**

By Alexander Goldin, PhD, Barbara Weinstein, PhD, and Nimrod Shiman

*The data show that each type of medical mask in this study essentially functioned as a low-pass acoustic filter for speech, attenuating the high frequencies (2000-7000 Hz) spoken by the wearer by 3 to 4 dB for a simple medical mask and close to 12 dB for the N95 masks. This means the speech quality*

### RECENT POSTS

HR Webinar: Hearing Healthcare Covid-19

The Musicians' Clinics of Canada Website

Duke Researchers Use Sound Waves to Transport Droplets on Rewritable Chip


CILcare Partners with Vertidialog Solutions

Interacoustics Academy to Host July 13-16

IHear Medical Launches Free @Hear Test

*degradation, in combination with room noise/reverberation and the absence of visual cues, renders speech close to unintelligible for many patients with hearing loss.*

The pervasive toll of hearing loss on clinician-patient communication has become magnified with advent of the coronavirus (COVID-19 virus). In fact, as is shown in **Figure 1**, the trajectory of hearing loss and the infection fatality by age distribution tell a very important story—one most of us may not have as of yet considered. In short, hospitals, emergency departments, and intensive care units (ICUs) are now flooded with older adults with the COVID-19 virus. Many of the people who have fallen victim to the virus have hearing loss, are unaccompanied by family members, are frail, have multiple chronic conditions and are likely without hearing assistance. Estimates of the infection fatality ratio in the oldest age group ( $\geq 80$  years) ranges from 8% to 36%, and sex-based data suggest that older men generally fare worse than elderly women in their ability to combat the virus.

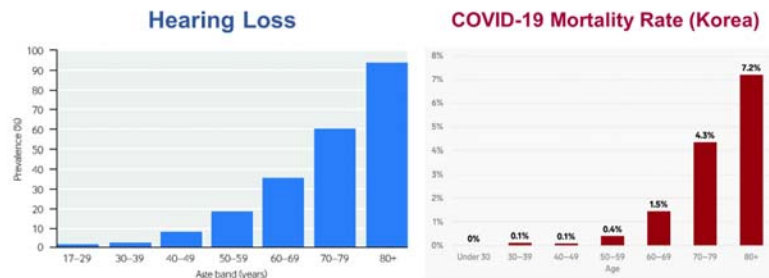

**Figure 1.** Hearing loss prevalence by age groups and mortality rates of persons by age groups in Korea with COVID-19 virus.

We are all aware that a cornerstone of patient safety and quality healthcare is effective communication which allows patients to participate fully in their care. For communication to be effective, it must take place in a manner appropriate to one's age, understanding, and communication abilities. Furthermore, medical information, when provided, must be **complete, accurate, timely, unambiguous, and understood** by the patient.<sup>1</sup> It is well established that when patients can hear and understand the healthcare provider, they are much more likely to be compliant and follow physician recommendations and have their medical needs optimally met.

Since social, psychological, and environmental factors also influence speech understanding especially in hospital settings, older adults are at a **great disadvantage** in what is now the new normal.<sup>2</sup> As is evident from **Table 1**, the hospital is a **very stressful listening situation**.<sup>3</sup> It is novel (N), unpredictable (U), threatening (T), and the skill set (S) to cope tends to be compromised. Add to this the fact that clinical care is now often being delivered in settings which are

#### FOLLOW US

 Facebook

 Twitter

 LinkedIn

 YouTube

extremely noisy and distracting—replete with alarms beeping and ongoing competing conversations between patients and providers.


Characteristic of a Stressful Listening Situation	Resilience Resource Being Challenged
Novelty	Coping
Unpredictability	Social Network Availability
Threat to Self	Emotional Awareness
Sense of Lack of Self Control	Self Efficacy

**Table 1.** *The Stressful Listening and Resilience Resources model, adapted from Pichora-Fuller (2017)<sup>3</sup>, can describe the challenges some patients will encounter during a hospital visit, particularly during the COVID-19 pandemic.*

Further, the struggle to understand speech for persons with hearing loss is now exacerbated by the masks which are a necessity to help minimize the spread of the virus. These masks pose two obvious problems for patients with hearing loss: 1) the patient cannot gain any cues from lipreading, and 2) the voice of the healthcare provider is attenuated and distorted. This brief study looks at the latter problem.


### How Does a Medical Mask Impact the Speech Signal?

We decided to quantify the potential challenge to a listener’s comprehension at the acoustic level posed by the variety of masks being worn by healthcare providers. The Alango Technologies team based in Israel measured voice samples in an anechoic chamber as a function of the type of mask being worn (**Figure 2**). A GRAS head and torso simulator was used to play white noise through the artificial mouth, and the output acoustic signal was measured by a microphone at 2 meters distance.

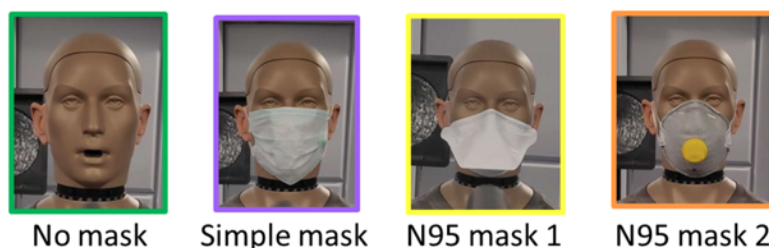

**Figure 2.** *Types of masks used in this study’s measurements.*

**Figure 3** summarizes the measurements. Only the differences between the green (no mask curve) and the others are important. Our findings were not unexpected but were unsettling. The data show that each mask essentially serves as a low-pass filter, attenuating the high frequencies (2000-7000 Hz) spoken by the wearer, with the decibel (dB) level of attenuation ranging from 3 to 4 dB for a simple medical mask and close to 12 dB for the N95 masks.

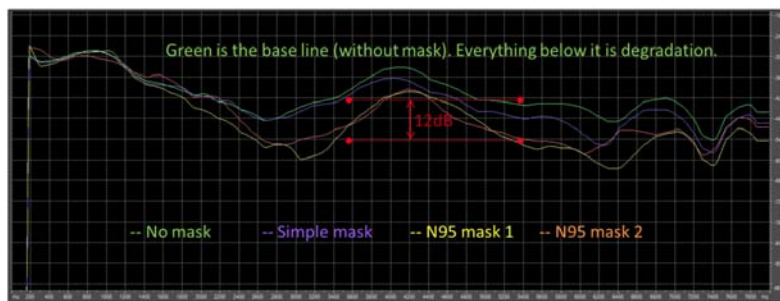

**Figure 3.** Output SPL level differences (in dB) by type of mask being worn.

In short, the speech quality degradation, in combination with room noise/reverberation and the absence of visual cues, renders speech close to unintelligible for many.

### Recommendations for Caregivers

We include in **Table 2** a list of universal communication precautions we have created, and encourage each of you to share with your patients and colleagues in healthcare settings. For those patients who do not have a hearing aid or assistive listening device, we also strongly recommend use of commercially available and affordable voice enhancement technologies (eg, Alango BeHear ACCESS or BeHear NOW, Williams Sound Pocket Talker, etc) or hearing apps on smartphones (eg, Jacoti ListenApp) to help optimize speech understanding.

Table 2. Tips for Talking to Patients in Hospital Settings
■ Speak slowly
■ Reduce noise and obtain the patient’s attention
■ Ask how patient prefers to communicate
■ Make sure hearing aid wearers are using them
■ Consider using a portable hearing amplifier
■ Rephrase remarks if not understood
■ Take turns while speaking to the patient
■ Do not shout, over-emphasize, or exaggerate your words
■ Do not speak to the patient while walking

**Table 2.** Tips for talking to patients in hospital settings

In the words of Helen Keller, it is our firm belief that “Alone we can do so little; together we can do so much.”

### References

1. Blustein J, Wallhagen MI, Weinstein BE, Chodosh J. Time to take hearing loss seriously. *Joint Commission Jour on Qual and Patient Safety.* 2020;46:53-58.
2. Peelle JE. Listening effort: How the cognitive consequences of acoustic challenge are reflected in brain and behavior. *Ear Hear.* 2018;39(2):204-214.

3. Pichora-Fuller MK. How social psychological factors may modulate auditory and cognitive functioning during listening. *Ear Hear.* 2016;37:92S-100S.

**About the authors: Alexander Goldin, PhD**, is the Founder and CEO of Alango Technologies in Tirat Carmel, Israel; **Barbara E. Weinstein, PhD**, is Professor & Founding Executive Officer of the Doctor of Audiology Program at City University of New York; and **Nimrod Shiman**, is an Audio Application Engineer at Alango Technologies.



**Alexander Goldin, PhD**



**Barbara Weinstein, PhD**

**CORRESPONDENCE** can be address to Dr Weinstein at: [bweinstein28@gmail.com](mailto:bweinstein28@gmail.com)

**Citation for this article:** Goldin A, Weinstein BE, Shiman N. How do medical masks degrade speech perception? *Hearing Review.* 2020;27(5):8-9.

SHARE:

RATE:


< PREVIOUS

NEXT >

[Audioscan Launches Updates for Verifit2 and Axiom Hearing Instrument Verification Systems](#)

[Sonic Expands Captivate Hearing Aid Family](#)

### RELATED POSTS

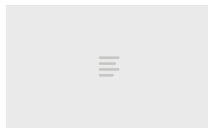


With over 200 million people worldwide suffering from hearing loss, hearing aids are a vital part of many people's lives. The hearing aid market is expected to reach \$10.5 billion by 2025.



**Richard...**



Hearing Charities of America is a national organization dedicated to providing hearing aids to those in need. We are currently seeking qualified individuals to join our team.

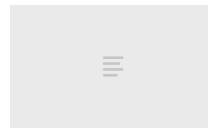

The University of Pittsburgh is seeking a qualified individual for a position in the Department of Hearing and Speech Sciences. The successful candidate will be responsible for teaching and research.